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
SMITH+NEPHEW Medical, Inc.
 Level 2 Zone 9, Daxing Park 15th Street, Suzhou Industrial Park, 215129 Suzhou, China
 UK: Smith+Nephew Ltd, 27 Gloucester Street, London, United Kingdom, WC1N 3BU
 USA: Smith+Nephew, Inc., 10000 North Central Expressway, Dallas, Texas 75243, USA
 Tel: +353 1 454 4444 Fax: +353 1 454 4445
 Tel: +1 817 381 8888 Fax: +1 817 381 8888
 Email: info@smith-nephew.com usa@smith-nephew.com
 Web: www.smith-nephew.com
 Tel: +41 76 551 11 11 Fax: +41 76 551 11 11
 Email: info@smith-nephew.ch
 Web: www.smith-nephew.ch

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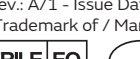
Winners Health Care, LTD.

10000 Highway 70, Unit 100, Langford, British Columbia, V1M 2B8, Canada
 Tel: +1 (250) 831-8888 Fax: +1 (250) 831-8888 Email: info@winnershealthcare.com
 Tel: +1 (800) 361-1111 Fax: +1 (800) 361-1111 Email: info@winnershealthcare.com

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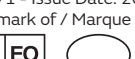the super absorbent dressing

Indication

per Absorbent Dressing is an absorbent pad designed for effective occlusive contact skin contact. This dressing is ideal for moderate to heavily exuding wounds. The unique super absorbent dressing can absorb up to 10 times its own weight in wound exudate. The dressing is designed to absorb and lock in wound exudate and prevent it from leaking out. The dressing is designed to be used on the skin and is not intended for use on the skin.

Instructions for Use

Application

1. Clean the wound area in accordance with normal procedures. Ensure the per-wound skin is dry.
2. The product is sterile before use. Do not use if the packaging is damaged. Select a dressing size to overlap the wound edges by at least 2-3cm.
3. Remove the film and then gently apply the dressing directly on the wound side with the adhesive side of the dressing onto the skin.
4. The dressing should be fixed with an adhesive tape. The dressing should be changed if the adhesive tape is damaged or if the adhesive tape is loose.
5. The dressing should be fixed off the skin the second day and remove the dressing.

Frequency of dressing change

The dressing should be changed every 7 days.

[illegible][illegible][illegible]